



Guangzhou Henco Electronic Technology Co., Ltd

HenCo[®] *for YOU*

MORE THAN PROFESSIONAL



Add: B16 ,4th floor Zhongxi Building No.306 Jiefangzhong Road
Yuexiu District Guangzhou China Zip:510030

Tel:0086-20-28326732

Mobil :0086-13760640202

Fax:0086-20-28326702

Web:www.hencobeauty.en.alibaba.com



HK-17 5 in 1 Tripolar RF Multipolar RF Vacuum Cavitation slimming machine

User Manual

RELEVANT INTRODUCTION

40 KHz ultrasonic principle

Powerful shock wave head emits frequency of 40 KHz sound waves into the human body. Which penetrates into human fat cells causes vibration impacts and movement in the fat cells and in turn can effectively burn calories, reduce water retention in fatty tissues and causes the fat cells to shrink. In addition, the strong acoustic vibration that causes cells to produce a strong vibration impact, cause fatty cells to bust instantaneous, thus reducing fatty cell numbers cells, allowing the removal of the fatty cells.

Radio Frequency Technology - RF

It is referred to as a high frequency AC change of electromagnetic waves. Change of more than 10,000 times a second called high-frequency current, while the RF is such a high frequency current. When RF waves penetrate the skin, it causes collagen fiber to heat up to 45°C-65°C, the collagen fibers immediately contract, so that the loose skin wrinkles are stretched by the stimulation of dermal collagen hyperplasia gradually promote cellular rejuvenation, so as to achieve the purpose of firming and reduce wrinkles.

Vacuum + bipolar RF theory

being integrated with the most advanced RF technology and radio frequency, the instrument may directly reach the deep-seated fat body and have the excellence of targeted positioning RF. In the fast-active state, fat cell tissue may generate friction heat, increase local temperature and remove excess fat and toxin from the body through sweat gland, enterohepatic circulation and lymph and ultimately achieve the effect of dissolving fat.

With controlled depth technology, inlaid diamond particles of different sizes may rub skin along its pattern, adjust suction strength by means of vacuum draw and rub with negative pressure strong force to directly explode thick fat. Thus, the effects are very obvious: Enhance blood circulation, accelerate the decomposition of fat, take away wastes, enhance the effect of removing cellulite, effectively improve the state of areolar tissue, remove stubborn fat and improve the sufficient absorption condition of diet products.

MAIN FUNCTIONS

- (1) Smooth fine wrinkles, shrink pores.
- (2) Make the skin moist.
- (3) Enhance lymphatic and blood circulation.
- (4) Get rid of facial redness.
- (5) Eliminate slow acne scars.
- (6) Promote collagen and cell activation.
- (7) Improve skin activity and toughness.
- (8) Increase the speed of metabolism, speed up the body to excrete waste and excessive water.
- (9) Reduce stretch marks.
- (10) Relax muscles, relieve muscle spasms, and relieve muscle pain.
- (11) To tighten muscles of the arms, legs, thighs, buttocks, lower back, abdominal muscles, re-shaping body contour.
- (12) Effectively improve the orange peel-like skin of the buttocks and thighs, while also helps in postpartum or after effect of liposuction in the abdomen area.

CAUTIONS

- (1) Before using the instrument, you should clean the skin, and apply the appropriate skin care products, to ensure that the instruments work properly.
- (2) Before using the instrument, remove metal jewelry removed so as not to weaken the function of the ion pulse.
- (3) People with sensitive skin will experience slight swelling, this is a normal phenomenon, and the swelling will soon subside.
- (4) 40K Fat Cavitation is recommended to be use once a week and two or three times per week.
- (5) Operation of the ultrasonic or 40K Fat Cavitation must have water and other skin care products; otherwise the instrument will have a certain loss.
- (6) 40K Fat Cavitation treatment for some customers may hear ringing in ears; this phenomenon is a normal phenomenon, because this band has been more than human hearing range.
- (7) This machine is 40 kHz shock wave fat Cavitation cannot be use on the head (neck). Not be used for the glottis and sexual organs.
- (8) After treatment, the user needs to drink much water to replenish moisture.

Do NOT SUITABLE TO THE PEOPLE BELOW

- (1) Individual with skin wound, open cuts
- (2) Epileptic suffers
- (3) Patients with heart with heart disease (arrhythmia), especially with a pacemaker
- (4) Cancer suffers
- (5) Patients with infectious diseases
- (6) Pregnant women
- (7) Major surgery patients, diabetic's sufferers
- (8) Metal implants patients.
- (9) While beginning the operation, ensure the handle is position at 90 degree (vertical)

to the skin; this is to avoid unease feeling by user


(10) To begin operation, first probe the head to the skin, and then switch on the main power switch. Switch the energy regulation of the mode of operation to a minimum, and then slowly increase the intensity.

INSTRUMENT SYSTEMS INSTRUCTION



Boot display and see Welcome screen, free to click on the screen to enter the function interface.


Main Interface

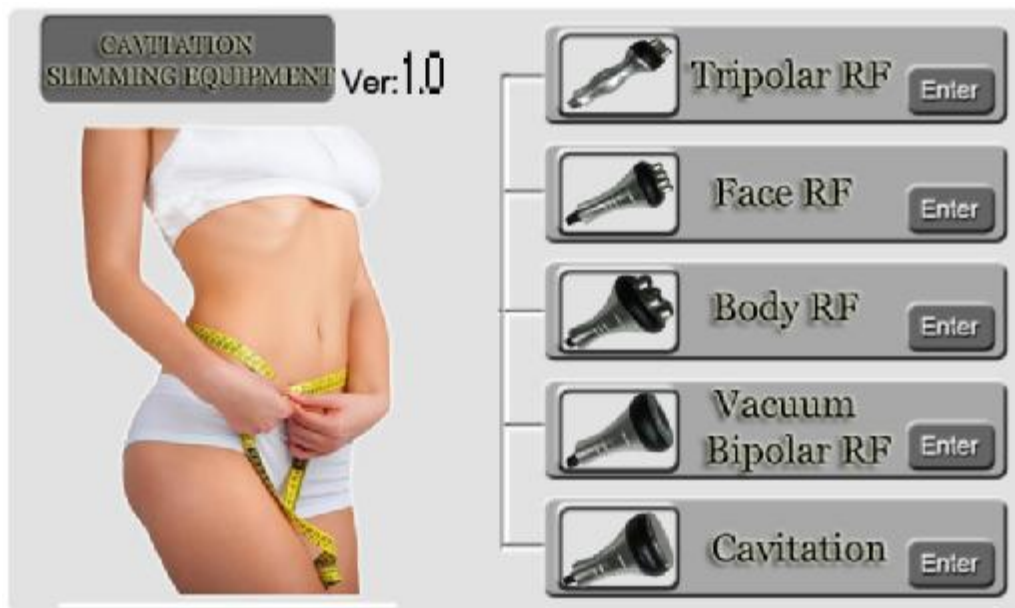

40K Fat Cavitation Function Interface




Body RF Function Interface

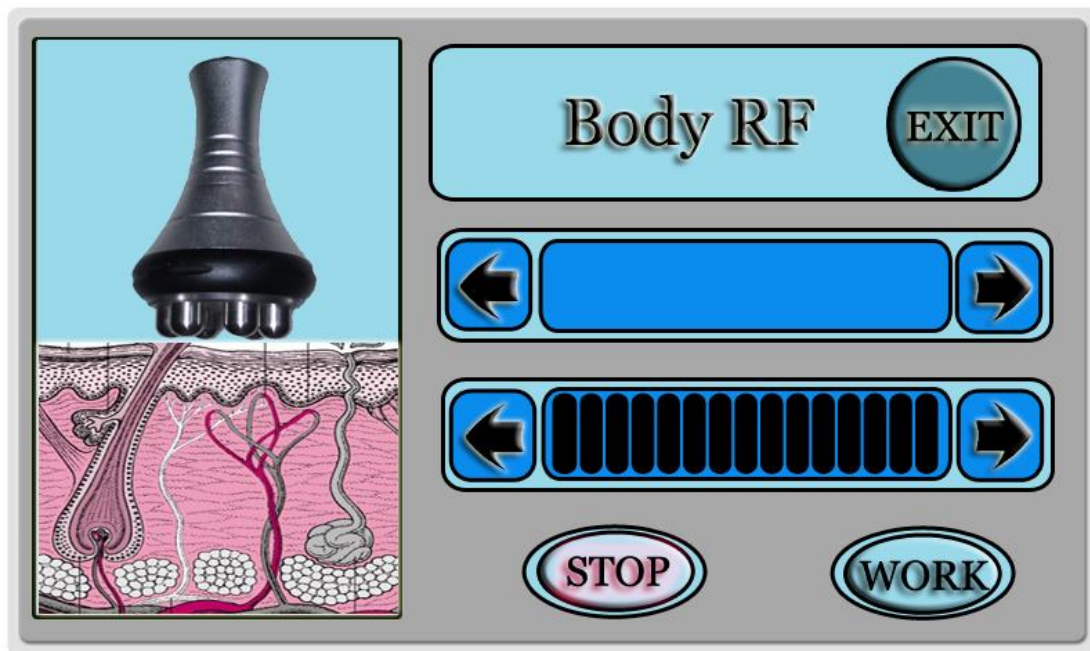


Face RF Function Interface



Tripolar RF Function Interface

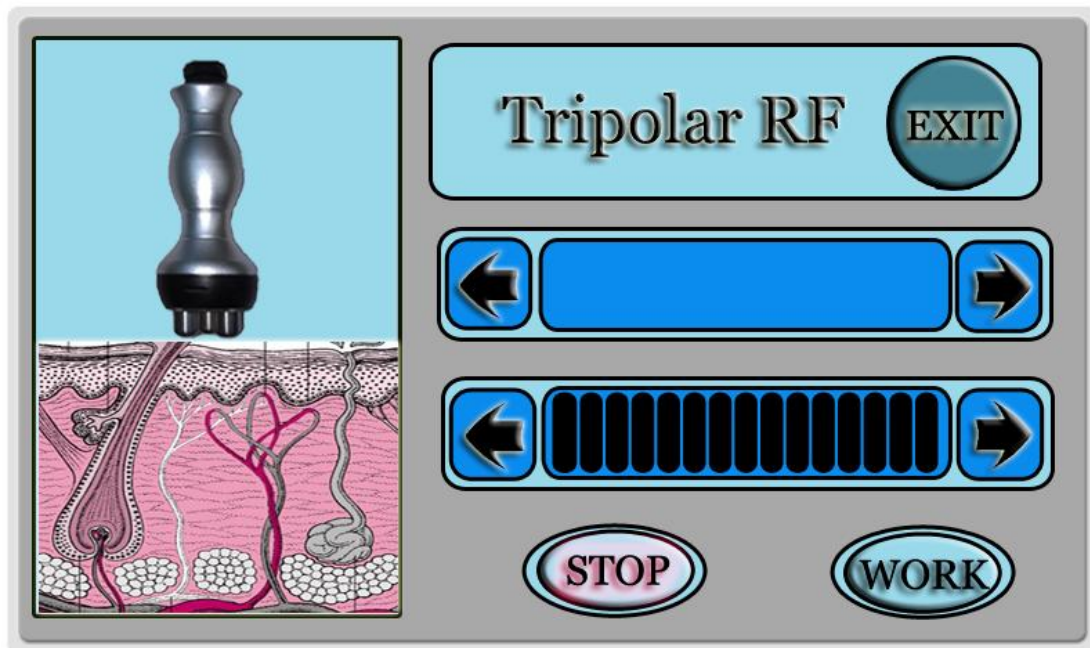

Vacuum+BipolarR Function Interface

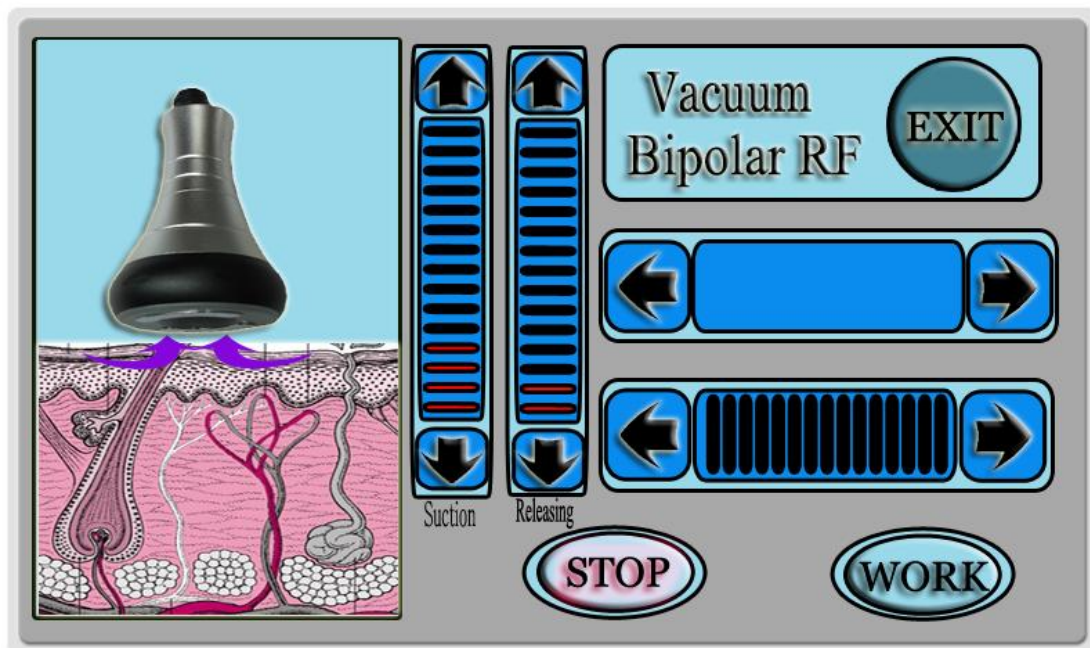

ACCESSORIES

Product Accessories



40KHz Cavitation



6-polar RF



4-polar RF



3-polar RF



Vacuum + RF

SUGGESTED OPERATING METHODS

Thigh

(1). From down to upper, pushing to the groin to dredge the lymph.




(2). From down to upper, circling push by anticlockwise can help decomposing fatness.



(3). Also do push by come-and-go to decompose fatness.



(4) Pull from the knee and muscle texture to upper, can improve the curve.

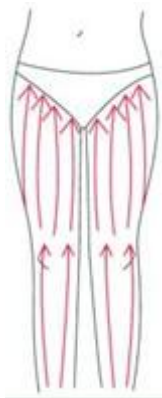


Buttock

(1) Along the lymph direction, pulling to the waist.



(2) From upper to down, by anticlockwise gesture, pull come-and-go motion.

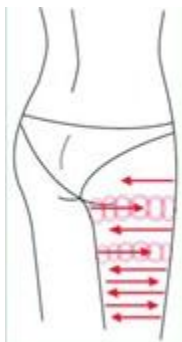


(3) In accordance with the muscle texture sliding compact buttocks skin. Shaping the hip circle slides back and forth. Perfect curve.

Back


(1) First by come-and-go, pushing the bladder nerve 2 to 3 times.



(2) Circling by anticlockwise to stimulate the underarm lymph node.

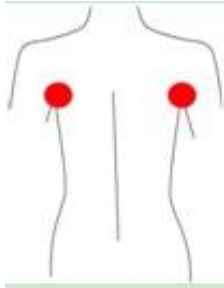

(3) On the back, by anticlockwise circling, can help to decompose the fatness.

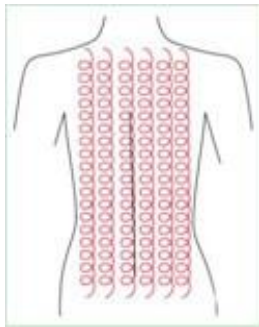

(4) Using lymph drainage gesture, pull the toxin to the lymph node, can help improving the back curve.



Arm

(1) From the inner elbow to the armpit, doing lymph drainage.



(2) Circling to stimulate the lymph node.

(3) By anticlockwise or come-and-go pushing, decomposing the fatness.




(4) From the elbow to the outer, doing tightening gesture.




Abdomen

(1) From small to large, by clockwise direction, circling around the navel can help peristalsis of the large intestine.

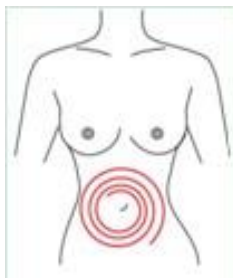

(2) From small to large, by anticlockwise direction, circling push down.

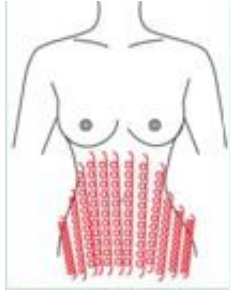

(3) From the belly and muscle texture, pull to the groin.



(4) Through lymph drainage gesture, taking the toxin to the place of groin.



SERVICE

If there is any question during the using process, please contact us.

Guarantee:

1. The warranty time is ONE year and lifetime maintenance.
2. Within warranty period, all exchanges or repair due to defects are free (Buyer is responsible for shipping charges).
3. If the product is out of the warranty time, our company is only responsible for maintaining, but the expenses fee and spare part fee should be charged.

WE ARE NOT RESPONSIBLE FOR THE FOLLOWING SITUATIONS:

1. The man-made damage or falls.
2. The problems are caused by disassembling, removing or remodeling the product by the user.
3. The damage caused by improper operation.

If any questions, please contact us promptly, and we will help you to solve the problem as soon as possible.

TECHNICAL SPECIFICATIONS

40K Cavitation:

Intensity: 50W/CM²

Max power consumption: ≤50W

Frequency: 40,000 Hz

Mode: firm/shape

Body radio frequency

Frequency: 1MHz

Radio frequency power: ≤40W

Face radio frequency

Frequency: 1MHz

Radio frequency power: ≤40W

Tripolar radio frequency

Frequency: 1 MHz

Radio frequency power: ≤40W

Vacuum + bipolar radio frequency

Frequency: 1 MHz

Radio frequency power: ≤40W

General information

Display: 7inch touch screen

Power: ≤80W

Operating voltage: 220V/50HZ; 110V/60HZ

Package size: 58.5cm*42.5cm*44.5 cm

G.W.:23kg